


# HPLC DATA SHEET

## Bisacodyl, Docusate Sodium Tablets (KP, Assay) (1)

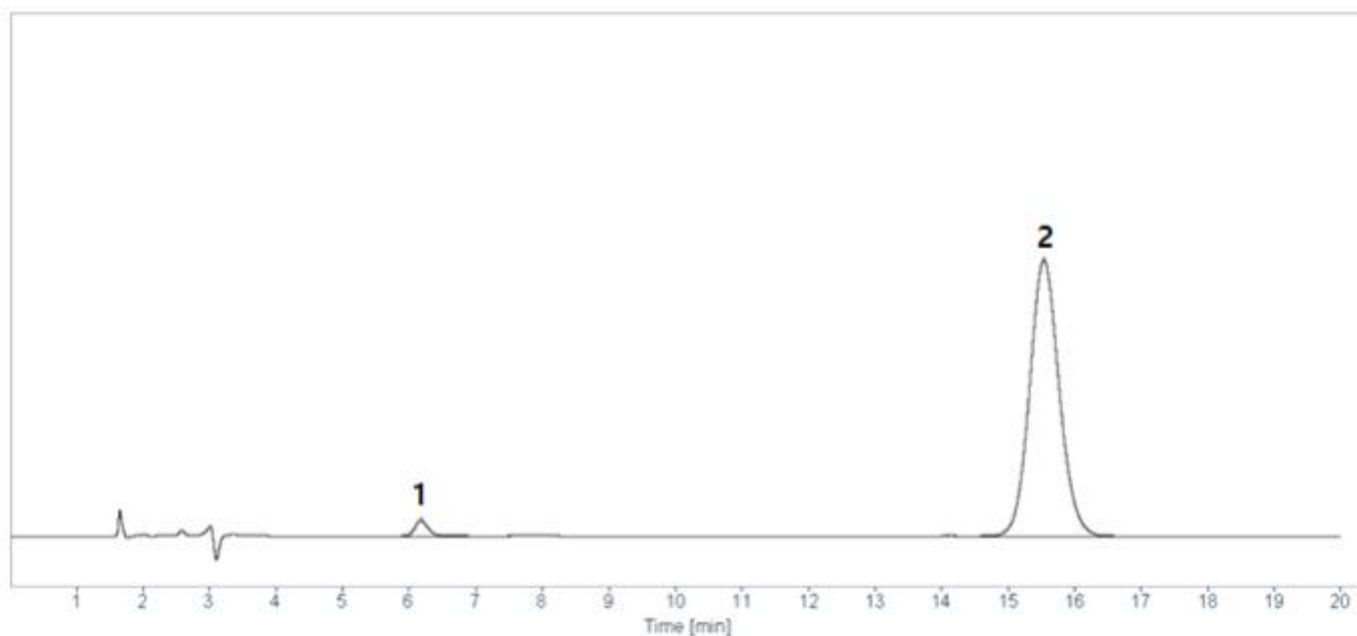

1. Docusate Sodium      2. Bisacodyl

### HPLC condition

Column            INNO Column C18 10um, 3.9mm I.d. × 300mm (P/N 010N03930)  
0.005 mol/L Tetrabutyl ammonium phosphate : Acetonitrile : Methanol = 26 : 23 : 1

Mobile Phase    – 0.005 mol/L Tetrabutylammonium phosphate : Dissolve 1.7 g of  
tetrabutylammonium phosphate in 1 L of water.

Flow rate        1.0 mL/min

Temperature     25 °C

Detector         UV (214 nm)

Injection Volume 10 uL