


HPLC DATA SHEET

Tocopherol Acetate (KP, Assay)

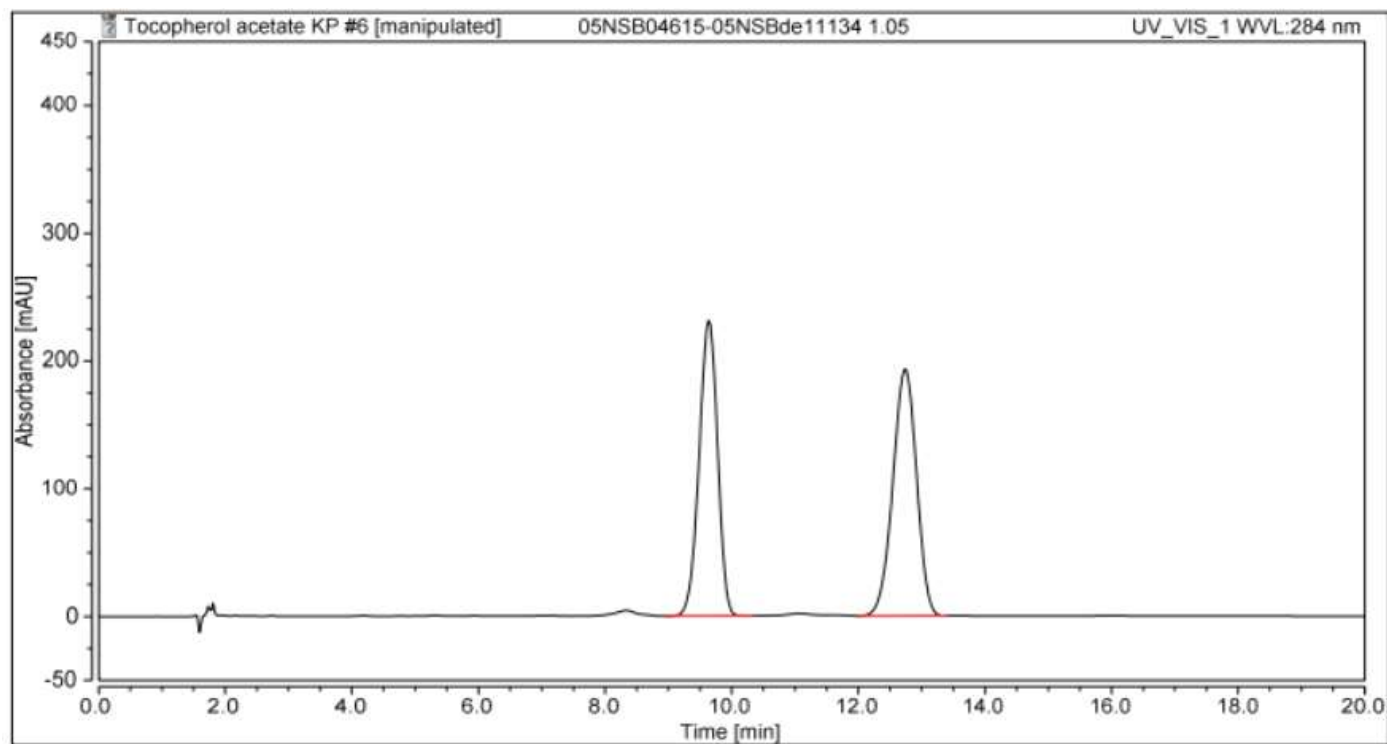

1. Tocopherol 2. Tocopherol Acetate

HPLC condition

Column	INNO Column C18-SB 5um, 4.6mm I.d. × 150mm (P/N = 05NSB04615)
Mobile Phase	Methanol : Water (49 : 1)
Flow rate	Adjust so that the retention time of Tocopherol Acetate is about 12 minutes (About 1.05 mL/min)
Temperature	35 °C
Detector	UV (284 nm)
Injection Volume	20 uL