


# HPLC DATA SHEET

Rebamipide (KP, 2-Amino-3- 2[(1H)-quinolinone-4-yl] propionate·2HCl)

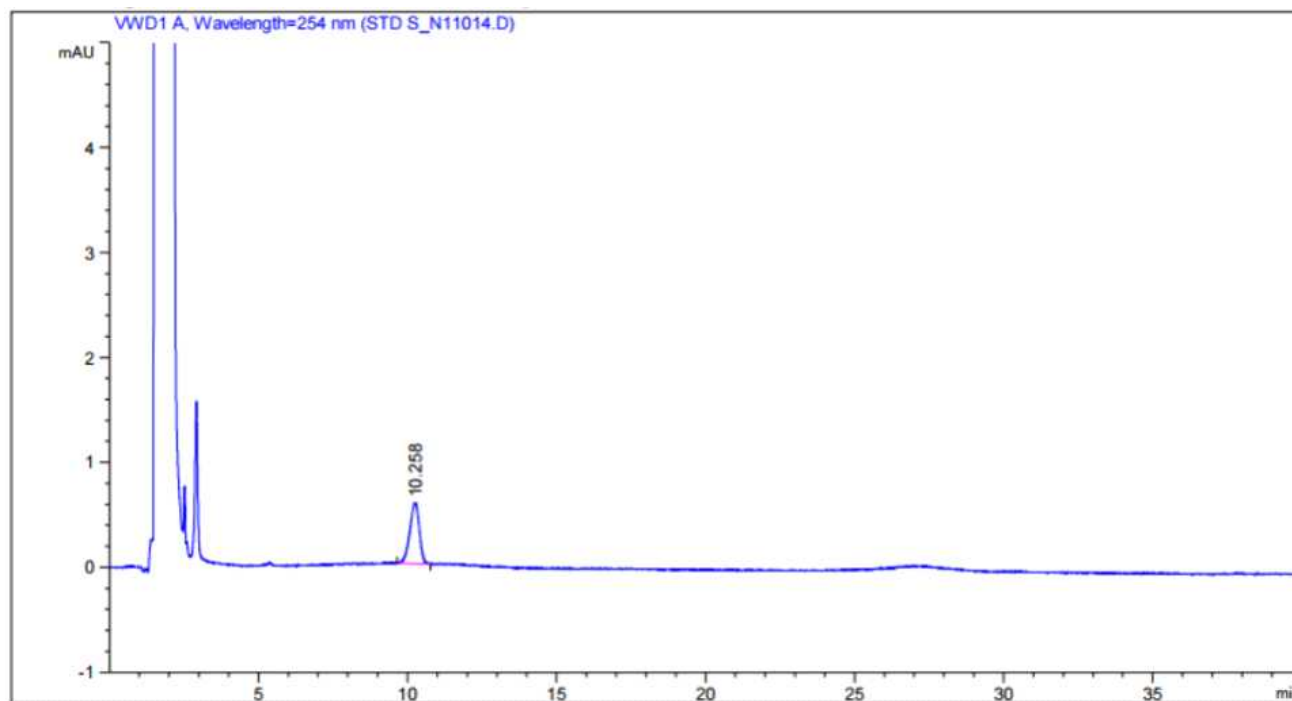

1. 2-Amino-3- 2[(1H)-quinolinone-4-yl] propionate·2HCl

## HPLC condition

|              |                                                                                                                                    |
|--------------|------------------------------------------------------------------------------------------------------------------------------------|
| Column       | INNO Column NH, 5um (4.0×150mm)                                                                                                    |
| Mobile Phase | Acetonitrile and Buffer (800 : 200)<br>* Buffer : 0.58g/L of sodium citrate dihydrate and 2.0g/L of anhydrous citric acid in water |
| Flow rate    | Adjust so that the retention time of is about 10 minutes.                                                                          |
| Temperature  | 25℃                                                                                                                                |
| Detector     | UV (250nm)                                                                                                                         |